

# beurer

## IPL hair removal



german|engineering

Our  
Beauty  
Our Beurer

Beautiful and well-groomed with the Beurer Beauty product range

Clinically tested technology for maximum safety at home.  
Safe, fast and reliable.

2) IPL 9000+ / IPL 10000+ dermatological study, 36 women, Institut proDERM Hamburg, study no.: 14.0148-11  
IPL 7500 / IPL Pure Skin Pro / IPL Velvet Skin Pro dermatological study, 20 women, Institut proDERM Hamburg, study no. 14.0309-11



## IPL long-lasting hair removal – developed by professionals


More and more women and men wish to have long-lasting smooth skin. Conventional methods such as shaving, epilation or waxing do not have the desired success. It does not take long for the unwanted hair to become visible again, and the applications – which are often painful – begin again. For this reason, those affected are now increasingly turning to long-lasting hair removal using intense pulsed light (IPL) devices. The light-based hair removal method is very fast and painless. Find out more about the topics of light-based hair removal, safety during application and get to know our fitting products:

## IPL technology for use at home

Beurer hair removal devices use the innovative “IPL” light-based technology, which is already used in professional beauty salons and by doctors for successful long-lasting hair removal. You can now use this tried and clinically tested technology safely and effectively at home. During use, our user-friendly IPL devices provide you with particular safety and protection of your skin. Save yourself the effort and the money for a trip to the beauty salon by using Beurer IPL devices.

## Our hair growth cycle

**To what extent does light-based hair removal impact on the hair growth cycle?** Each hair on the body goes through three phases of the hair growth cycle: the *anagen phase* (growth phase), the *catagen phase* (transitional phase) and the *telogen phase* (resting phase).

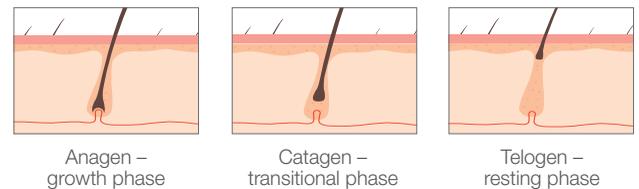


The length of a full hair growth cycle varies from person to person and depends on the body area. As a rule, a cycle is normally 18 to 24 months. Only hairs that are in the anagen phase respond to treatment with the IPL technology. Therefore, it is important to note that multiple treatments are needed to treat all hairs. At least one complete hair growth cycle is needed to achieve long-lasting hair removal results.


## Safe, fast and long-term – clinically tested technology

The Beurer IPL devices feature a safety system with integrated UV protection and 2-in-1 skin type and skin-contact sensor.


The skin-contact sensor protects your eyes by ensuring that light pulses are only emitted when the light area is in full contact with your skin. The integrated UV filter also protects your skin from harmful rays.

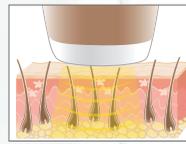

Before application

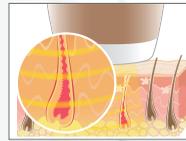

During a pulse

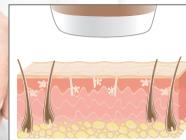

After application



Light opening with integrated UV filter

Integrated 2-in-1 skin type and skin contact sensor for the utmost safety

With useful “Auto Flash” mode for convenient application

With skin type recognition – automatically sets the right energy level



## The principle behind light-based hair removal

How can light remove hair? Light-based hair removal is founded on the principle of selective photothermolysis, which uses optical energy to deactivate hair growth. This effect is achieved by selective absorption of the light energy and its conversion into heat.

Melanin is the pigment in the hair shaft that is responsible for absorbing the light and thereby generating the heat that deactivates the hair follicle deep under the skin and prevents hair growth. The more melanin there is in the hair (i.e. the darker the hair is), the more light can be absorbed and the more effective the hair removal is. The IPL devices remove the hair using light pulses which work over a wavelength range of 475–1200 nm.



Convenient IPL technology for use at home



Our IPL devices can be used on (naturally) dark blond, brown and black hair and on light to medium skin tones. They are not effective on very light or grey hair. The long-lasting hair removal<sup>3)</sup> with IPL is also not suitable for naturally dark skin tones or particularly tanned skin.

To ensure protection, our devices are equipped with a skin type sensor that measures skin pigmentation on the surface to be treated before each light pulse. For your safety, the device then only triggers if the skin pigmentation is suitable.

Skin colour						
Hair colour						
	X	X	X	X	X	X
	✓	✓	✓	✓	✓	X

## Application period

**Body treatment:** Just 3-4 treatments can achieve a hair reduction of up to 50 %. For body treatments, you can achieve the best results with 4 treatments carried out every 2 weeks, followed by 3 treatments carried out every 4 weeks. The treatment can then be repeated as required.



## Application areas

The IPL treatment is intended for use by adults for the removal of unwanted body hair.

The ideal body areas are legs, armpits, bikini line, chest, abdomen and face below the cheekbones for women. The device is not intended for removing hair on the anus, genitals or nipples.



Face



Bikini line



Armpits



Legs



Arms



Chest and back


<sup>3)</sup> Long-lasting results over several months

“beurer MyIPL” app – for additional support before, during and after treatment

Effective, quick and convenient – these are the features of the Beurer IPL hair removal devices. The new “beurer MyIPL” app helps you to create a personalised treatment plan for long-lasting hair removal.

#### What the app offers:


- Safe professional use at home
- Energy level recommendation for the treatment
- Personalised treatment plan through the precise coordination of personal hair and skin characteristics (hair removal regions)
- Calendar view and reminder function
- Safety notes on correct usage of the IPL device
- Organisation of treatment dates throughout the entire period of use

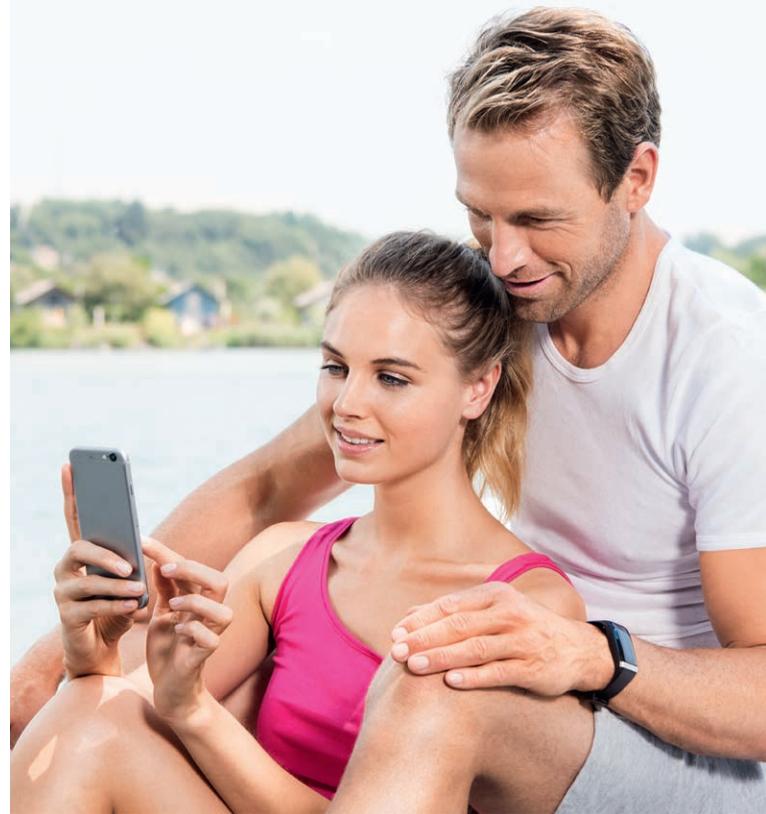


#### Create your personal treatment plan in just a few steps:

1. Download the free “beurer MyIPL” app
2. Select the Beurer IPL hair removal device and the area of the body to be treated
3. Determine the hair and skin colour type
4. Automatic creation of a personal treatment plan

“beurer MyIPL” app – modern design, intuitive operation and helpful tutorial

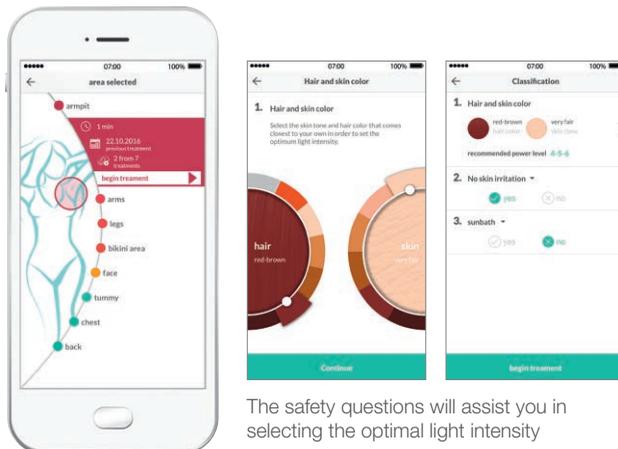

The safety questions will assist you in selecting the optimal light intensity

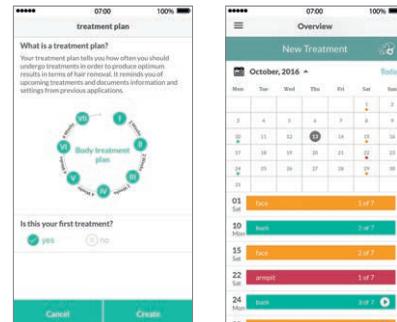

To achieve the best possible result, the treatment plan informs and reminds you of the correct application intervals



# Which IPL hair removal device is suitable for me?

The appeal of our devices for long-lasting hair removal lies in their simple and safe operation for use at home. You benefit from our innovative professional technology, which is not only reliable and safe, but also clinically tested and dermatologically approved for skin tolerance<sup>2)</sup>.

- Long-lasting hair removal with clinically tested professional technology
- State-of-the-art light technology used by dermatologists for professional hair removal
- Maximum safety with 2-in-1 skin type and skin contact sensor
- For the face<sup>3)</sup>, arms, legs, armpits, bikini line and the back, chest and abdomen
- Including integrated UV filter to protect the skin



## IPL 7500 SatinSkin Pro

- Compact device – ideal for smaller body areas
- Up to 150,000 light pulses
- 3.1 cm<sup>2</sup> light surface for precise application
- Flash-glide function for uninterrupted application
- 3 energy levels
- Mains operation



## IPL Pure Skin Pro

- Compact device – perfect for treating small and hard-to-reach areas of the body
- Up to 200,000 light pulses
- 3.1 cm<sup>2</sup> light surface for precise application
- “Auto Flash” mode for convenient application
- 3 energy levels
- Mains operation



## IPL Velvet Skin Pro

- Battery-powered device – incl. precision attachment for face and bikini line
- Up to 300,000 light pulses
- 4.5 cm<sup>2</sup> light surface for fast application
- 2.0 cm<sup>2</sup> precision attachment
- “Auto Flash” mode for convenient application
- 6 energy levels



## IPL 9000+ SalonPro System

- Hair removal device – for the fastest treatment on the market
- Up to 100,000 light pulses
- Replacement cartridge with up to 100,000 light pulses available
- XXL – large 7 cm<sup>2</sup> light surface for particularly fast application
- Flash-glide function for uninterrupted application
- 6 energy levels
- Mains operation



## IPL 10000+ SalonPro System

- Hair removal device – for the fastest treatment on the market
- Incl. precision attachment for the face
- Lifetime flashes – up to 250,000 light pulses
- XXL – large 7 cm<sup>2</sup> light surface for particularly fast application
- With auto sensor for individual skin tone detection
- Flash-glide function for uninterrupted application
- 6 energy levels
- Mains operation

2) IPL 9000+ / IPL 10000+ dermatological study, 36 women, Institut proDERM Hamburg, study no.: 14.0148-11  
IPL 7500 / IPL Pure Skin Pro / IPL Velvet Skin Pro dermatological study, 20 women, Institut proDERM Hamburg, study no. 14.0309-11  
3) On the face only for women, below the cheek line



## Skin tolerance dermatologically approved<sup>2)</sup>

Long-lasting hair removal for use at home – the dermatologically approved IPL technology guarantees safety and skin tolerance from the first application.



<sup>2)</sup> IPL 9000+ / IPL 10000+ dermatological study, 36 women, Institut proDERM Hamburg, study no.: 14.0148-11  
IPL 7500 / IPL Pure Skin Pro / IPL Velvet Skin Pro dermatological study, 20 women, Institut proDERM Hamburg, study no. 14.0309-11

Beurer GmbH. P.O. Box 1427. 89004 Ulm, Germany. Soeflinger Str. 218. 89077 Ulm, Germany  
Tel. +49 (0) 731/3989-0. Fax. +49 (0) 731/3989-139. [export@beurer.de](mailto:export@beurer.de)  
[www.beurer.com](http://www.beurer.com) [www.beurer-healthguide.com](http://www.beurer-healthguide.com)

