



HealthManager
Pro

beurer
connect

beurer HealthManager Pro

beurer
CONNECT



Easily obtain a holistic picture of your health and keep an eye on it at all times




Medical device app that meets the requirements of the relevant European directives



Download the "beurer HealthManager Pro" app now free of charge!

HealthManager Pro

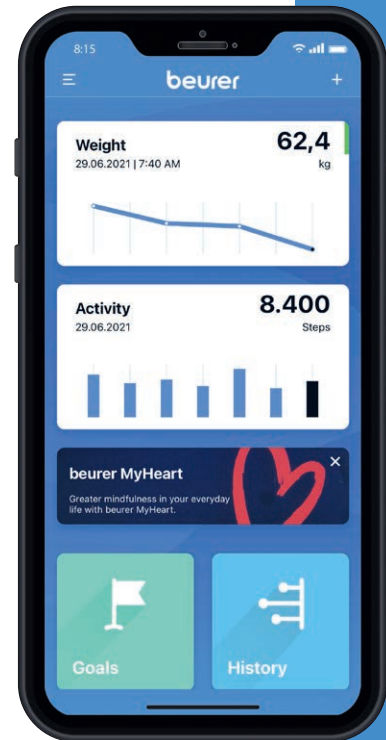

All-in-one solution: App can be combined with over 30 Beurer Connect products from the following categories: activity, blood pressure, blood glucose, ECG, weight, pulse oximetry, temperature and drinking

Medical device app that meets the requirements of the relevant European directives

Modern design and intuitive operation

Transfer app data: Transfer data from the old "beurer HealthManager" app to the new "beurer HealthManager Pro" app

"beurer MyHeart": Program to help on the way to a healthier lifestyle (additional service subject to a charge)

"beurer MyCardio Pro": Easy analysis of ECG measurements from home (additional service subject to a charge)

Accessible: Large areas to click, easy-to-read fonts and high contrasts

Optimised, guided device integration process

Measured values presented clearly and in detail in tables and graphics

Function to export data to a PDF or CSV file to easily exchange information with your doctor

Measured values are classified based on reference values or personal goals

Managing medication in the "Medication cabinet" area


Comment function for noting down certain information, e.g. extreme values

Data relating to weight, blood pressure, activity, sleep, temperature and pulse oximeter can be shared with Apple Health

Data relating to activity and sleep can be shared with Google Fit

Multi-user function: Several users can use one app

App is available in 22 languages: German, English, French, Spanish, Italian, Greek, Dutch, Czech, Polish, Portuguese, Swedish, Finnish, Slovak, Hungarian, Danish, Romanian, Slovenian, Russian, Traditional Chinese, Simplified Chinese, Turkish, Bulgarian

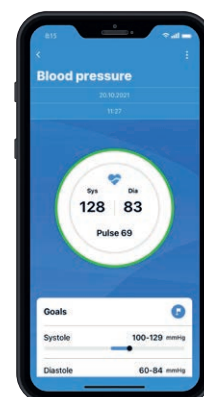

HealthManager-Pro_2022-06-12_01_DB_BEU_EN Subject to errors and changes



Additional services in “beurer HealthManager Pro”


beurer MyHeart

The holistic “beurer MyHeart” program helps to integrate a healthier lifestyle into everyday life.

The four elements: healthy recipes, exercise, knowledge and daily inspiration guide through the journey for 30 days towards a healthy future.

Additional service subject to a charge.

Can be booked for 30 days or for an unlimited period.

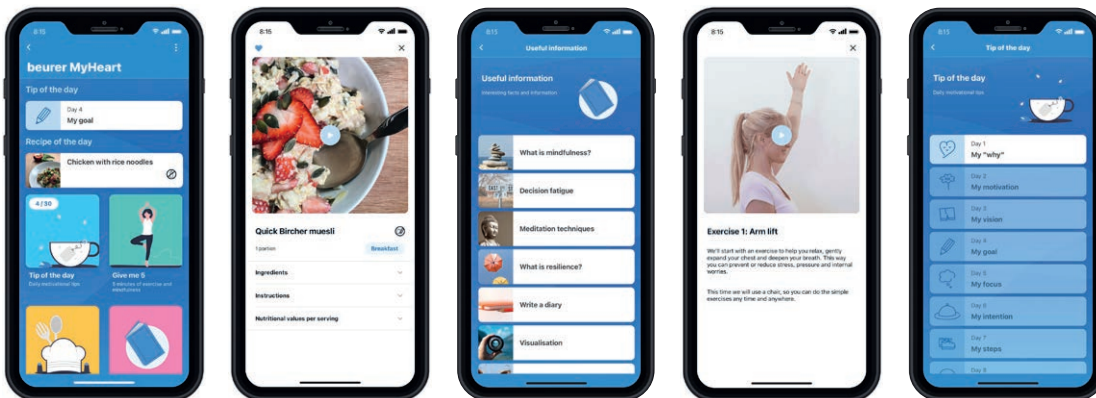


beurer MyCardio Pro

Taking ECG measurements at home and obtaining an analysis has never been easier. Simply taking the ECG measurement in familiar surroundings and transferring the measured values to the “beurer HealthManager Pro” app.

Using the “beurer MyCardio Pro” service, a detailed analysis of the ECG measurements is immediately created, along with a professional report for sending to the doctor.

Additional service subject to a charge.

Can be booked in packages of 1 report, 10 reports & 100 reports.

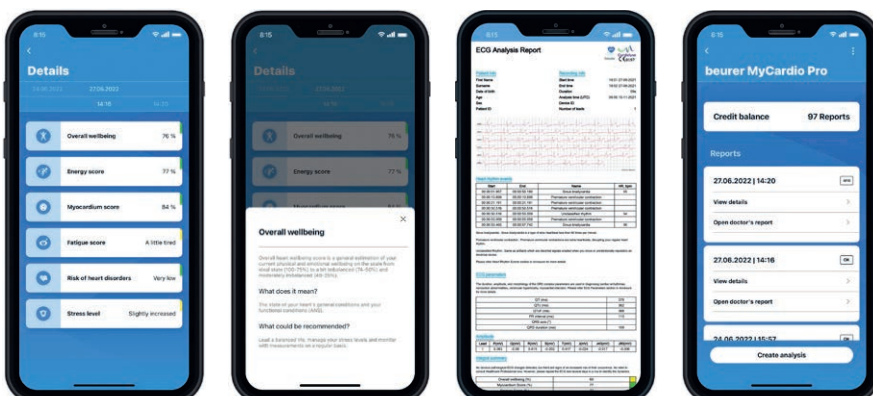




Transfer app data easily and without hassle

App data can easily be transferred from the old "beurer HealthManager" app to the new "beurer HealthManager Pro" app.


1

- Open the settings in the **"beurer HealthManager"** app
- Select the menu item **"Transfer app data"**

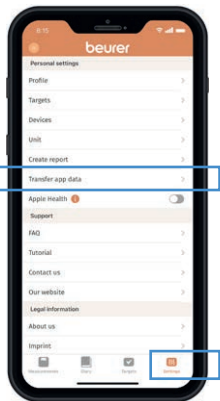


2

- Tap the **"beurer HealthManager Pro"** app logo
- Install the app

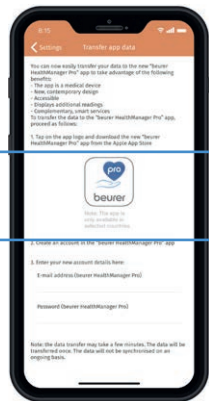


3

- Open **"beurer HealthManager Pro"** and register again
- **Save your account online**

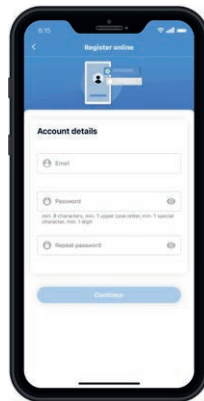

4

- Return to **"beurer HealthManager"** Under the **"Transfer app data"** menu item, enter the account data of **"beurer HealthManager Pro"** that you previously created and transfer the data
→ Done!

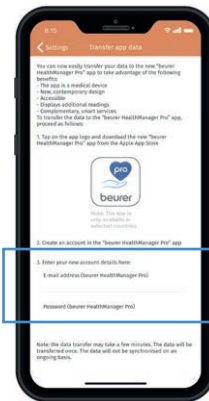

works with



Synchronisation with
iOS ≥ 14.0
Android™ ≥ 8.0



List of supported devices